

## **Drug & Medicine Dispensing**

124	Dispensing froneys
129	CLIQ® & RFID Solutions
130	Patient Administration Trolleys
132	Drawer Trolleys
133	Tray Trolleys
134	Blister Pack (MDS) Trolleys
136	Nomad® Cassette Trolleys
138	Biodose® & MultiMeds Trolleys
139	Medinoxx Trolleys
140	Laptop Holders
142	Security Clamp
143	Distribution Trolleys

Diamanaina Trallava



CTDISPM1/10/3S2D/HSB

Item	Description	Weight
High Security Bolt Lock		
CTDISPM1/10/3S2D/HSB	3 x 100mm, 2 x 150mm Trays, High Security Bolt Lock	35.5kg
CTDISPM1/10/4S1D/HSB	4 x 100mm, 1 x 150mm Trays, High Security Bolt Lock	35.2kg
CTDISPM1/10/6S/HSB	6 x 100mm Trays, High Security Bolt Lock	36.3kg

#### **Caretray Trolley - Dispensing**

#### **Features**

- Trolleys suitable for the storage & dispensing of drugs & medicines
- Suitable for use with Electronic Prescribing & Medicines Administration (EPMA)
   & Electronic Patient Records (EPR)
- · Quick & easy disassembly, assists thorough & deep cleaning
- Fully welded mild steel frame, powder coated in a durable epoxy paint (requires no on-site assembly)
- Recyclable polypropylene one-piece top & side panels, smooth surfaces & radiused edges assists cleaning (compatible with steam cleaning)
  - Top Panel Raised edges on three sides provide defined grab points to act as push handles, whilst an integral rim at the front contains spillages (max load 30kg)
  - Side Panels Incorporate runners with positive stops preventing trays from being accidentally withdrawn
- · Carerails fitted to either side support a wide range of accessories
- Solid grade laminate door & back panel providing security for the contents of the travs
- Two locking options...
  - · High Security Bolt Lock Supplied with two keys
  - Electronic Push Button Lock Programable four-digit code
- Trolley pre-configured with the following accessories...
  - · Waste Bin & Bracket, 8 Litre
- Translucent, impact-resistant trays & dividers, smooth surfaces & radiused edges assists cleaning (compatible with steam cleaning / high temperature machine washing)
- Trays feature an integral handle with a plastic clip-on label cover (coloured labels supplied - white, blue, green, yellow & red)
- 100 / 150mm deep trays, can be subdivided in both length & width using removable dividers (max load 5kg)
- 2 x long & 2 x short dividers supplied with each tray
- Polyurethane ring buffers offer protection to both the trolley & fabric of the building, smooth surfaces assist cleaning
- · 100mm polymer castors, 2 off braking
- When not in use, trolleys storing medicines should be secured to the wall using a Security Clamp (see locking & security accessories)

#### Dimensions (w x d x h)

#### **Standard Colours**

- Overall Dimensions 550 x 550 x 993mm
- Plastic Top 455 x 455 (internal) at height 968mm
- Trays
  - Shallow Tray 355 x 405 x 100mm (internal)
  - Deep Tray 355 x 405 x 150mm (internal)

· Paint - Spearmint

#### **Options / Accessories**

- Locking & Security page 60
- Overbridge & Tilt Bins page 66
  Carerail Mounted page 68

- Laptop Holders page 74
   Castors & Ring Buffers page 76
   Trays, Shelves, Dividers & Labels page 77



Item	Description	Weight
Electronic Push Button L	ock	
CTDISPM1/10/3S2D/EPB	3 x 100mm, 2 x 150mm Trays, Electronic Push Button Lock	35.5kg
CTDISPM1/10/4S1D/EPB	4 x 100mm, 1 x 150mm Trays, Electronic Push Button Lock	35.2kg
CTDISPM1/10/6S/FPB	6 x 100mm Trays Electronic Push Button Lock	36 3ka



CTDISPM1/10/3S2D/EPB

www.bristolmaid.com



Item	Description	Weight	
High Security Bo	olt Lock		
MT100	Small, High Security Bolt Lock	26.0kg	
MT110	Medium, High Security Bolt Lock	35.0kg	
MT120	Large, High Security Bolt Lock	42.5kg	
Electronic Push Button Lock			
MT100CL	Small, Electronic Push Button Lock	26.0kg	
MT110CL	Medium, Electronic Push Button Lock	35.0kg	
MT120CL	Large, Electronic Push Button Lock	42.5kg	

#### **Dispensing Trolleys - Wooden, Standard Capacity**

#### **Features**

- Trollevs suitable for the storage & dispensing of drugs & medicines Lockable hinged lids...
- Upper Assisted using a gas strut, providing document support

  - Lower Provides a work area, at height 980mm (max load 4.0kg)
- Two locking options...
  - High Security Bolt Lock Supplied with two keys
  - Electronic Push Button Programmable four-digit code
- Hinged work flap provides a work area at height 930mm (max load 2.0kg)
- Lower shelf
- Five tier plastic rack (395mm)
- · Two slide out storage / waste bins
- 100mm pressed steel castors
- · When not in use, trolley should be secured to the building see optional MT/1/6F

#### Dimensions (w x d x h)

- Overall...
  - Small 750 x 490 x 1130mm
  - Medium 890 x 490 x 1130mm
  - Large 1030 x 525 x 1130mm
- Internal...
  - Small 530 x 427 x 263mm sloping to 112mm at front
  - Medium 683 x 427 x 263mm sloping to 112mm at front
  - Large 805 x 461 x 263mm sloping to 112mm at front

#### **Standard Colours**

- Frame Grey
- · Shelves / Dividers White
- · Cabinet Spearmint
- Shelf / Flap Grev

#### **Options / Accessories**

- MT/1/6F Security Wall Clamp
- · ZP005 Ring Buffers

sales@bristolmaid.com +44 (0)1258 484455 www.bristolmaid.com



sales@bristolmaid.com

+44 (0)1258 484455

Item	Description	Weight
High Secur	ity Bolt Lock	
MT220	Medium, High Security Bolt Lock	42.0kg
MT230	Large, High Security Bolt Lock	50.0kg
Electronic	Push Button Lock	
MT220CL	Medium, Electronic Push Button Lock	42.0kg
MT230CL	Large, Electronic Push Button Lock	50.0kg

#### **Dispensing Trolleys - Wooden, High Capacity**

#### **Features**

- Trolleys suitable for the storage & dispensing of drugs & medicines
- Lockable hinged lids...
  - · Upper Assisted using a gas strut, providing document support
  - Lower Provides a work area, at height 921mm (max load 4.0kg)
- · Two locking options...
  - High Security Bolt Lock Supplied with two keys
  - Electronic Push Button Programmable four-digit code
- Hinged work flap provides a work area at height 850mm (max load 2.0kg)
- Lower shelf
- Five tier plastic rack (395mm)
- · Two slide out storage / waste bins
- 125mm pressed steel castors
- When not in use, trolley should be secured to the building see optional MT/1/6F

#### Dimensions (w x d x h)

- · Overall...
  - Medium 760 x 560 x 1290mm
  - Large 905 x 560 x 1290mm
- · Internal...
  - Medium 530 x 461 x 513mm sloping to 137mm at the front
  - Large 683 x 461 x 513mm sloping to 137mm at the front

#### **Standard Colours**

- · Frame Grey
- Shelves / Dividers White
- Cabinet Spearmint
- Shelf / Flap Grey

#### **Options / Accessories**

• MT/1/6F - Security Wall Clamp

www.bristolmaid.com



Item	Description	Weight
MT130	Small, High Security Bolt Lock	22.0kg
MT140	Medium, High Security Bolt Lock	25.0kg
MT150	Large, High Security Bolt Lock	31.0kg

#### **Dispensing Trolleys - Metal, Standard Capacity**

#### **Features**

- Trolleys suitable for the storage & dispensing of drugs & medicines
- Lockable hinged lid provides document support
- Fitted with a high security bolt lock (supplied with two keys)
- Hinged writing flap provides a work area at height 907mm (max capacity 2.0kg)
- Lower shelf
- Five tier plastic rack (395mm)
- · Two slide out storage / waste bins
- 100mm pressed steel castors
- When not in use, trolley should be secured to the building see optional MT/1/6F


#### Dimensions (w x d x h)

- · Overall...
  - Small 750 x 459 x 1112mm
  - Medium 891 x 459 x 1112mm
  - Large 1030 x 492 x 1112mm
- Internal...
  - Small 489 x 457 x 292mm sloping to 127mm at the front
  - Medium 706 x 457 x 292mm sloping to 127mm at the front
  - Large 832 x 490 x 292mm sloping to 127mm at the front

#### **Standard Colours**

- · Frame Grey
- · Cabinet White
- · Shelf / Flap Spearmint

- MT/1/6F Security Wall Clamp
- · ZP005 Ring Buffers

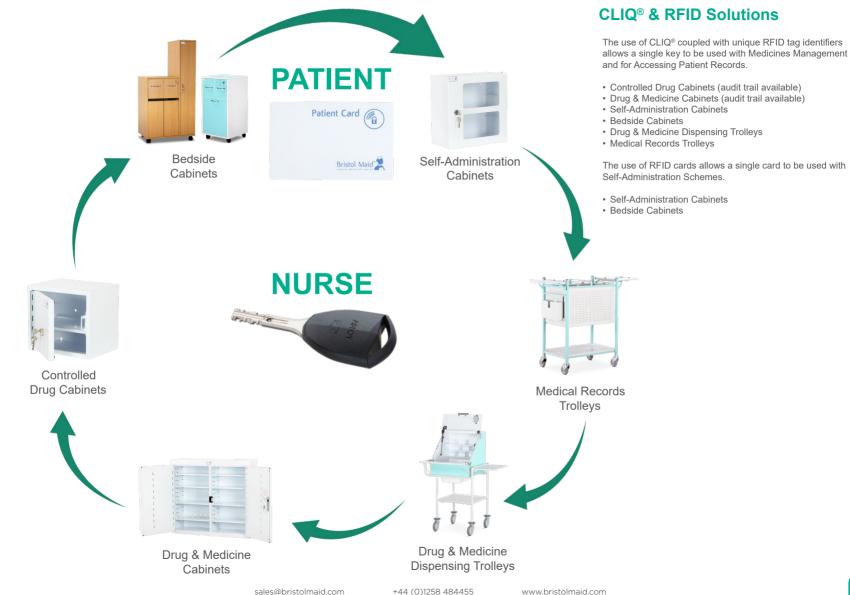



#### **Single Door**

High Security Bolt Lock	Electronic Push Button Lock	Description	Weight
PTS/HSB/PA24	PTS/EPB/PA24	Patient Administration, 24 Compartments	51.4kg

#### **Pharmacy Trolleys - Patient Administration**

#### **Features**

- Trollevs suitable for the storage & dispensing of drugs & medicines using individual compartments dedicated for each patient
- · Two locking options...
  - · High Security Bolt Lock Supplied with two keys
  - Electronic Push Button Programmable four-digit code
- Supplied with 4 / 8 trays subdivided into equal compartments...
  - 24 Compartment 4 travs subdivided into 6 equal compartments, each compartment - 190mm x 132mm x 220mm (5.5 Litres). Label holders provided for each compartment
  - 32 Compartment 8 trays subdivided into 4 equal compartments, each compartment - 192mm x 143mm x 220mm (6.1 Litres). Label holders provided for each compartment
- · Upper work surface
- Push handles incorporate integral buffers to offer protection to both the trolley & fabric of the building
- · Plastic debris bin & bracket
- · 125mm polymer castors
- When not in use, trolley should be secured to the building see optional MT/1/6F

#### Dimensions (w x d x h)

#### **Standard Colours**

- Single Door 627 x 579 x 1165mm
- · Paint Grey
- Double Door 820 x 579 x 1165mm

- · Work Surface Spearmint

- MT/1/6F Security Wall Clamp
- PT/LHR/S Laptop Holder Removable, Swing Arm
- PT/LHNR/S Laptop Holder Non-Removable, Swing Arm



PTS/EPB/PA24



Trays tilt to assist access to contents



PTS/HSB/PA24



PTD/HSB/PA32

#### **Double Door**

High Security Bolt Lock	Electronic Push Button Lock	Description	Weight
PTD/HSB/PA32	PTD/EPB/PA32	Patient Administration, 32 Compartments	68.8kg



Item	Description	Weight
PTS/HSB/DR6	High Security Bolt Lock, Six Drawers	53.0kg
PTS/EPB/DR6	Electronic Push Button Lock, Six Drawers	53.0kg

#### **Pharmacy Trolleys - Drawers**

#### **Features**

- Trolleys suitable for the storage & dispensing of drugs & medicines using drawers
- Two locking options...
  - High Security Bolt Lock Supplied with two keys
  - Electronic Push Button Programmable four-digit code
- Upper work surface
- Push handles incorporate integral buffers to offer protection to both the trolley & fabric of the building
- 3 adjustable door shelves
- · Plastic debris bin & bracket
- 125mm polymer castors
- When not in use, trolley should be secured to the building see optional MT/1/6F

#### Dimensions (w x d x h)

- Overall Dimensions 511 x 542 x 1072mm
- Drawer 335 x 295 x 113mm (internal)

#### **Standard Colours**

- · Paint Grey
- · Work Surface Spearmint

- MT/1/6F Security Wall Clamp
- PT/LHR/S Laptop Holder Removable, Swing Arm
- PT/LHNR/S Laptop Holder Non-Removable, Swing Arm



PTS/EPB/DR6



High Security Bolt Lock	Electronic Push Button Lock	Description	Weight
Single Door			
PTS/HSB/LA16	PTS/EPB/LA16	Sixteen LA Trays	39.0kg
Double Door			
PTD/HSB/A18C12	PTD/EPB/A18C12	Eighteen A & Twelve C Trays	57.0kg
PTD/HSB/LP24	PTD/EPB/LP24	Twenty Four LP Trays	57.0kg

#### **Pharmacy Trolleys - Trays**

#### **Features**

- Trolleys suitable for the storage & dispensing of drugs & medicines using individual patient trays
- · Two locking options...
  - High Security Bolt Lock Supplied with two keys
  - · Electronic Push Button Programmable four-digit code
- · Trays supported on adjustable shelves
- Tray Sizes (w x d x h)...
  - LA 86 x 268 x 57mm (aluminium)
  - A 85 x 181 x 54mm (plastic)
  - C 129 x 181 x 54mm (plastic)
  - LP 115 x 266 x 95mm (plastic)
- · Upper work surface
- Push handles incorporate integral buffers to offer protection to both the trolley & fabric of the building
- · Adjustable door shelves...
  - Single Door 3 door shelves
  - Double Door 6 door shelves
- · Plastic debris bin & bracket
- 125mm polymer castors
- When not in use, trolley should be secured to the building see optional MT/1/6F

#### Dimensions (w x d x h)

- Single Door 511 x 542 x 1072mm
- Double Door 711 x 568 x 1072mm

#### **Standard Colours**

- Paint Grey
- · Work Surface Spearmint

#### **Options / Accessories**

- MT/1/6F Security Wall Clamp
- PT/LHR/S Laptop Holder Removable, Swing Arm
- PT/LHNR/S Laptop Holder Non-Removable, Swing Arm

sales@bristolmaid.com +44 (0)1258 484455 www.bristolmaid.com

# PTS/HSB/BP6 PTS/HSB/DR/BP6 (frames not included) PTS/HSB/BP4 (frames not included)

#### **Single Door**

High Security Bolt Lock	Electronic Push Button Lock	Description	Weight
PTS/HSB/BP4	PTS/EPB/BP4	Blister Pack, Four Frames	31.0kg
PTS/HSB/BP6	PTS/EPB/BP6	Blister Pack, Six Frames	38.0kg
PTS/HSB/DR/BP6	PTS/EPB/DR/BP6	Drawer, Blister Pack, Six Frames	39.0kg

#### **Pharmacy Trolleys - Blister Pack (MDS)**

#### **Features**

- Trollevs suitable for the storage & dispensing of drugs & medicines using the blister packed monitored dosage system (MDS)
- · Two locking options...
  - High Security Bolt Lock Supplied with two keys
  - Electronic Push Button Programmable four-digit code
- Frames supported on either a fixed panel mounted to the inside rear of the trolley or on three adjustable slide out shelves (frames not supplied)
- Upper work surface
- Push handles incorporate integral buffers to offer protection to both the trolley & fabric of the building
- · Adjustable door shelves...
  - Single Door 3 door shelves
  - Double Door 6 door shelves
- Fold over MDS dispensing surface
- Plastic debris bin & bracket
- 125mm polymer castors
- · When not in use, trolley should be secured to the building see optional MT/1/6F

#### Dimensions (w x d x h)

- Single Door 511 x 542 x 1072mm
- Double Door 711 x 568 x 1072mm

#### **Standard Colours**

- · Paint Grey
- · Work Surface Spearmint

- MT/1/6F Security Wall Clamp
- PT/LHR/S Laptop Holder Removable, Swing Arm
- PT/LHNR/S Laptop Holder Non-Removable, Swing Arm



PTS/EPB/BP4



PTS/HSB/DR/BP6

PTD/HSB/DR/BP9





PTD/HSB/DR/BP9

PTD/HSB/BP6 (frames not included)

#### **Double Door**

High Security Bolt Lock	Electronic Push Button Lock	Description	Weight
PTD/HSB/BP6	PTD/EPB/BP6	Blister Pack, Six Frames	38.0kg
PTD/HSB/BP9	PTD/EPB/BP9	Blister Pack, Nine Frames	42.0kg
PTD/HSB/DR/BP9	PTD/EPB/DR/BP9	Drawer, Blister Pack, Nine Frames	46.0kg

Tab prevents the shelf from being removed

135



#### **Single Door**

High Security Bolt Lock	Electronic Push Button Lock	Description	Weight
PTS/HSB/NC30	PTS/EPB/NC30	Nomad®, Thirty Cassettes	38.0kg

#### **Pharmacy Trolleys - Nomad® Cassettes**

#### Features

- Trolleys suitable for the storage & dispensing of drugs & medicines using the Nomad® cassette monitored dosage system
- Two locking options...
  - High Security Bolt Lock Supplied with two keys
  - Electronic Push Button Programmable four-digit code
- Cassettes supported on three adjustable shelves (cassettes not supplied)
- Upper work surface
- Push handles incorporate integral buffers to offer protection to both the trolley & fabric of the building
- Adjustable door shelves...
  - Single Door 3 door shelves
  - Double Door 6 door shelves
- · Plastic debris bin & bracket
- · 125mm polymer castors
- When not in use, trolley should be secured to the building see optional MT/1/6F

#### Dimensions (w x d x h)

- Single Door 511 x 542 x 1072mm
- Double Door 711 x 568 x 1072mm

#### **Standard Colours**

- · Paint Grey
- · Work Surface Spearmint

#### **Options / Accessories**

- MT/1/6F Security Wall Clamp
- PT/LHR/S Laptop Holder Removable, Swing Arm
- PT/LHNR/S Laptop Holder Non-Removable, Swing Arm



PTS/EPB/NC30

sales@bristolmaid.com +44 (0)1258 484455 www.bristolmaid.com





#### **Double Door**

High Security Bolt Lock	Electronic Push Button Lock	Description	Weight
PTD/HSB/NC56	PTD/EPB/NC56	Nomad®. Fifty Six Cassettes	48.0ka

www.bristolmaid.com

#### Pharmacy Trolleys - Biodose® & MultiMeds

# PTD/HSB/BD69 (trays not included)

High Security Bolt Lock	Electronic Push Button Lock	Description	Weight
Single Door			
PTS/HSB/BD46	PTS/EPB/BD46	Biodose® & MultiMeds, 46 Trays	37.3kg
Double Door			
PTD/HSB/BD69	PTD/EPB/BD69	Biodose® & MultiMeds, 69 Trays	52.0kg

PTS/EPB/BD46 (travs not included)

#### **Features**

- Trolleys suitable for the storage & dispensing of drugs & medicines using the Biodose® & MultiMeds monitored dosage system
- · Two locking options...
  - High Security Bolt Lock Supplied with two keys
  - Electronic Push Button Programmable four-digit code
- Trays supported on fixed runners (trays not supplied)
- Upper work surface
- Push handles incorporate integral buffers to offer protection to both the trolley & fabric of the building
- · Adjustable door shelves...
  - Single Door 3 door shelves
- Double Door 6 door shelves
- · Plastic debris bin & bracket
- 125mm polymer castors
- When not in use, trolley should be secured to the building see optional MT/1/6F

#### Dimensions (w x d x h)

- Single Door 581 x 568 x 1072mm
- Double Door 812 x 568 x 1072mm

#### **Standard Colours**

- · Paint Grev
- · Work Surface Spearmint

- MT/1/6F Security Wall Clamp
- PT/LHR/S Laptop Holder Removable, Swing Arm
- PT/LHNR/S Laptop Holder Non-Removable, Swing Arm



High Security Bolt Lock	Electronic Push Button Lock	Description	Weight
Single Door			
PTS/HSB/MX40	PTS/EPB/MX40	Medinoxx, 40 Trays	41.4kg
Double Door			
PTD/HSB/MX60	PTD/EPB/MX60	Medinoxx, 60 Trays	55.2kg

#### **Pharmacy Trolleys - Medinoxx**

#### **Features**

- Trolleys suitable for the storage & dispensing of drugs & medicines using the Medinoxx monitored dosage system (compatible with 7 x 4 medi-trays)
- Two locking options...
  - · High Security Bolt Lock Supplied with two keys
  - Electronic Push Button Programmable four-digit code
- Trays supported on fixed runners (trays not supplied)
- · Upper work surface
- Push handles incorporate integral buffers to protect the fabric of the buildina
- Adjustable door shelves
  - · Single Door 3 door shelves
  - · Double Door 6 door shelves
- · Plastic debris bin & bracket
- 125mm polymer castors offer reduced rolling resistance making the trolley easier to manoeuvre, lowers noise levels, whilst smooth surfaces assist cleaning
- · When not in use, trolley should be secured to the building - see optional MT/1/6F

#### Dimensions (w x d x h)

- Single Door 581 x 568 x 1072mm
- Double Door 812 x 568 x 1072mm

#### **Standard Colours**

- · Paint Grey
- · Work Surface Spearmint

- MT/1/6F Security Wall Clamp
- PT/LHR/S Laptop Holder Removable, Swing Arm
- PT/LHNR/S Laptop Holder Non-Removable, Swing Arm



PTS/EPB/MX40



#### **Laptop Holder - Removable, Swing Arm**

#### **Features**

- · Three section swing arm offering height, rotation & angle adjustment
- Compatible with the PT range of single & double door pharmacy trolleys
- Fitted to an upright bolted to the rear left-hand corner of the trolley (upright & holder supplied for customer to fit)
- Laptop supported in a cradle, fitted with adjustable wing arms
- Laptop can be secured using a security clamp (supplied)\*
- Adjustable to accept laptops 280 435 x 22mm (w x h)
- · Maximum capacity 2.0kg
- · Design allows easy access to laptop ports

\*Not compatible with Dell laptops using the N17 wedge lock system

#### Standard Colours

· Paint - Grey

Item	Description	Weight
PT/LHR/S	Removable, Swing Arm	4.7ka

#### Laptop Holder - Non-Removable, Swing Arm

#### **Features**


- Three section swing arm offering height, rotation & angle adjustment
- Compatible with the PT range of single & double door pharmacy trollevs
- Fitted to an upright bolted to the rear left-hand corner of the trolley (upright & holder supplied for customer to fit)
- Laptop supported in a cradle fitted with a steel mounting bracket, fixed using either nyloc or wing nuts
- Laptop can be secured using a security clamp (supplied)\*
- Adjustable to accept laptops 280 435 x 22mm (w x h)
- · Maximum capacity 2.0kg
- · Design allows easy access to laptop ports

\*Not compatible with Dell laptops using the N17 wedge lock system

#### **Standard Colours**

· Paint - Grev

Item	Description	Weight
PT/LHNR/S	Non-Removable, Swing Arm	5.0kg

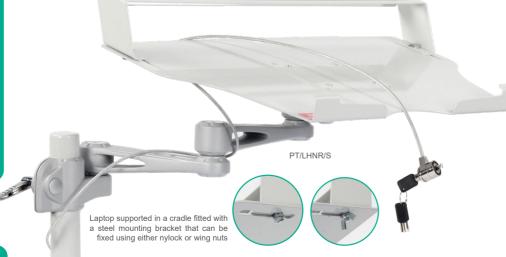

sales@bristolmaid.com +44 (0)1258 484455 www.bristolmaid.com






#### **Fitting Your Security Clamp**



Securely fasten one end of the clamp to the wall

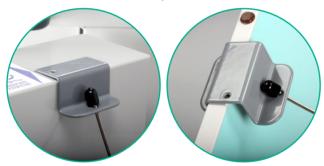

The 'hook shaped' end is trapped by the door when closed & locked



#### **Security Clamp**

#### **Features**

- · Allows the trolley to be secured to the 'fabric of the building'
- · Requires no additional keys, utilises the trolleys own locking system
- One end is secured to the wall, the 'hook shaped' end is trapped by the door/lid of the trolley when closed & locked
- · Recommended wall fitting height approx 1000mm



Important Note - Security Clamps

It is a requirement that when the trolley is not in use it should be secured to the 'fabric of the building'. Not only does this prevent casual theft, i.e. trolleys being pushed to an isolated area & broken into, but also encourages good nursing practice, giving the trolley a "permanent home".

"Safe & Secure Handling of Medicines - Royal Pharmaceutical Society, GB-March 2005"

Item	Description	Cable Length	Weight
MT/1/6F	Security Clamp	950mm	0.75kg



### ItemDescriptionWeightSA005Distribution Trolley - 15 Drawers49.0kgSA025Distribution Trolley - 18 Drawers56.0kg

#### **Distribution Trolleys - Lockable Drawers**

#### Features

- Trolleys suitable for the distribution of drugs & medicines from pharmacy to the ward & also with self-administration / patient own drug schemes
- Provides secure drug delivery using drawers which when locked are securely fastened to the trolley
- Drawers are fitted with a CAM lock (supplied with two keys)
- Max load per drawer 3kg
- Supplied as standard with locks to differ, allowing each ward to have its own key & be responsible for removing the drawer from the trolley (remaining key should be kept at the pharmacy)
- Master keys available (additional cost)
- Upper work surface
- Push handles at either end
- Document pocket fitted to end of trolley
- Low level corner bumpers offer protection to both the trolley & fabric of the building
- · 125mm pressed steel castors, 2 off braking

#### Dimensions (w x d x h)

- SA005 1026 x 523 x 1002mm
- SA025 1180 x 523 x 1002mm
- Drawer 145 x 423 x 230mm (Internal)

#### **Standard Colours**

- Paint Grev
- · Work Surface Spearmint

#### **Options / Accessories**

· Master Keys